


Electric microcrystallograph – SC632



USER MANUAL

Parts Name

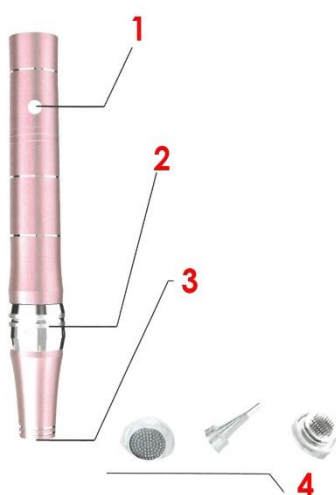

1. Adjustment of speed and ON/OFF; 2. Adjustment of needle length; 3. Needle connect; 4. Needle cartridge


How to Operate

1. Usage method

- 1) Clean skin: daily make-up removal, cleansing, can use toner for secondary cleaning.
- 2) Take out the nanowafer and install it on the top of the microcrystallizer.
- 3) Apply proper amount of liquid or essence to your face.
- 4) Press the power switch key for 3 seconds, start the microcrystallograph, and contact the skin vertically. Each contact point is recommended to stay for 35 seconds.

5) After the operation, please turn off the power supply and remove the nano chip, apply the common skin care products, and recommend the use of essence rice. Make a sterile repair mask before applying skin care products.

6) Nano chip can not be used repeatedly, so as to avoid two times of pollution.

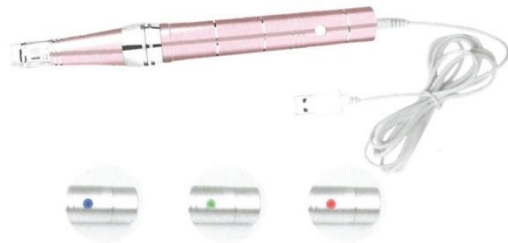

Instrument operation

1) Long press start/shut down, default first blue light is on, click the button to switch gear cyclically after start, second green light is on, third red light is on.


2) Charging display red light, full charge display green light.

How to Operate

Preparation

(1) Take out needle cartridge in the sterilized packing.

(2) Fitting round slots 1 into 2 of machine, push needle cartridge horizontally up to the end.

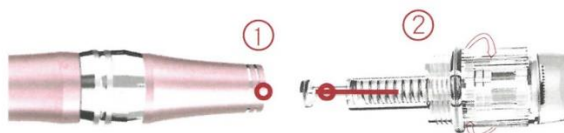

(3) Lock the cartridge by turning it to the counter clockwise

(4) Adjustment of needle length should be done after turn on the machine

(5) If you want to remove the needle cartridge, turn it to the clockwise and push out horizontally, fitting round slots (slot 1 and slot 2) together.

Adjustment of needle length

You can adjust the length of needle by turning the adjustment ring clockwise or counter clockwise. Turning it clockwise, length of it will be longer. Turning it counter clockwise, length of it will be shorter.



1. Counter clockwise: Length of it will be shorter.
2. Clockwise: Length of it will be longer.

Manipulation techniques

Point method: When the microcrystallograph is installed, the nanowafers are vertically touched on the skin surface, and a slight depression can be observed visually, with a gentle and dense upper and lower edge. Adjust the strength according to the contact position

Sliding method: When the microcrystallograph is installed, the microcrystalline wafer is vertically touched on the skin surface, and a slight quadruple-pit sensation can be observed visually. While sliding in a straight line, the nanowafer does not leave the skin surface. It is recommended to stay at the same part for 3-5 seconds and no more than 2 times.

Circle method: When the microcrystallograph is installed, the nano-wafer is perpendicular to the skin surface. It can visually detect the sense of retraction. The nano-wafer is circled at a uniform speed and does not leave the skin surface. It is recommended that the nano-wafer stay at the same place for 3-5 seconds and no more than 2 times.

Matters needing attention

Notes for nano-electric microcrystallograph:

- 1) Please do not disassemble or refit chargers and microcrystals. It is strictly forbidden to use chargers when the power supply is damaged first. Otherwise, it will lead to electric shock, fire and charger damage.
- 2) It is strictly forbidden to contact charger with wet hands. If charger touches other liquids in water and fire, it should cut off the power supply immediately.
- 3) No short-circuit charger, no mechanical vibration or shaking charger, no direct sunlight, no bathroom and other high-temperature places to use charger, the equipment is not waterproof, no rain or damp.
- 4) Please do not use chargers near TV, radio and other electrical appliances.
- 5) Pay attention to the charger when plugging in and pulling the power cord will damage the wire. When charging, the microcrystallizer and charger must be placed in the ventilation place with the temperature of - 10 C + 55 C. Using the unauthorized charger can cause danger, and also violate the warranty terms stipulated by the manufacturer.
- 6) Please keep the charger out of the reach of children to avoid danger.

Notes for nanoneedles

- 1) If the disposable packing is incomplete and leaky, please do not use it to avoid skin contamination.
- 2) Before operation, the operating parts and hands should be clean and clean.
- 3) For the first time, it is recommended to do a skin test with iskin care products on the inside of the wrist. Do not use it if itching and other allergic . phenomena occur.

4) The function of the needle to exert moderate force, not to press intentionally, nor to lift the wrist intentionally until the wafer can not even touch the skin.

5) About 10-20 minutes after the use of nano-wafers, some people will have skin redness, which is a normal phenomenon and will disappear within 24 hours without treatment.

6) Except for the needle head, no contact with water or other liquids is allowed in other parts of the instrument in order to avoid damaging the machine.

7) The same needle can only be used by one person for one time. It is strictly forbidden to share one wafer with others. It is forbidden to use the same wafer many times after disinfection. Otherwise, it may cause skin allergies, red feet, spots, inflammation and other adverse consequences.

8) When preserving nanoneedles, discard used wafers and store unopened needles and microcrystallographs in a well-ventilated environment to avoid dampness, heavy pressure, and shock, and stay away from children.

9) In the process of use, the nano-wafer should move continuously and stay at the same part for 3-5 seconds and no more than 2 times.

Suggestions for the use of Dry and sensitive skin:

Before using all kinds of efficacy products, first three import products (hyaluronic acid replenishment mask and after skin does not need water, then effective products will be imported.

Before introducing effective products, small-scale introduction should be made. Skin is not red, swollen, itchy and other symptoms, and then used in the face.

Please use the wash-free repairing and rehydrating products in time after importation.

Suggestions for skin with spot problems:

Freckle essence combined with replenishment products, such as whitening and freckle essence, is used.

Skin with acne problems: It can be imported with acne products, acne skin itself has inflammation, it is recommended not to use the whole face back and forth when importing, only in the place of acne, so as to avoid infection of other parts of the face.

Print Import Note: It can be imported with acne removing products and freeze-dried powder. First use freeze-dried powder to import the pockmarks. When importing, avoid infection by passing the operation of the pockmarks, finally, introduce the anti-acne products vertically on the Sockmarks.

Specific manipulation of acne pit: After cleaning the face, first use the instrument to point on the pit, then apply freeze-dried powder on the pit, and then use the instrument to point until the pit is red or slightly painful. Two times a week, daily morning and evening smear frozen thousand powder, generally continue to use for more than three months, the effect is more obvious.