



Massage warm palace belt  
(Chargeable)

manual




EMS microelectronics technology: through the micro-current AI intelligent control chip to stimulate the movement of fat cells to achieve the effect of weight loss and weight loss

2. EMS pulse massage instrument function introduction



Eight modes, 22 dynamics

3. Precautions for use

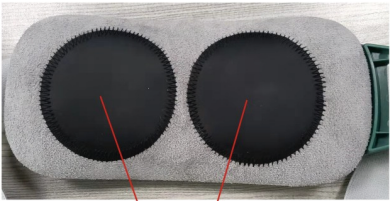

When using the **massage mode**, it can be used normally after being moistened with water (EMS) in the future (**it must be close to the belly, not across the belly**)

4.Product parameter

product name	sp1801
Input voltage	USB/DC5.1V
Power consumption	0.5W
Rechargeable lithium battery	3.7V/1800MAH
Battery Life	2 years
Input Current	Maximum 15MA
Output current	1.2MA
Working status	Working current is 5MA
Frequency output to the patch	1-100HZ
Product Features	Fever vibration pulse massage
Controller unit	4 gears
Motor vibration mode	4 kinds
Automatic light machine function	It will automatically shut down within 30 seconds
operating hours	It will automatically shut down after 15 minutes
charging	The red light goes out after fully charged
Charging port	USB charging port

5. Instructions for use

1. Vibration massage needs to be attached to the belly to operate normally, and the heating can be heated through the clothes.
2. It will automatically shut down after 15 minutes of smart timing
3. **Intelligent protection** can not work when charging, preventing the pulse magnetic field from chaos during charging, reducing risk, safety and guarantee
4. **Intelligent power off**, it will automatically power off after fully charged, and the battery will not respond when it is charged
5. There are a total of **four strengths and five temperatures** to choose from
6. Weight loss requires time-effective fat cell memory cycle is 30 days, weight loss basis is greater than 30 days, stick to 30 days slimming effect will be better
7. **Do not exceed 1 hour** a day, you need to add water in time

6.Safety Precautions

- one. Do not use this product with the following hospital electronic equipment colleagues, otherwise it may cause an accident or cause the medical electronic equipment to malfunction and cause serious damage to the body.
1. Implant wet medical electronic devices such as pacemakers
  - 2 Artificial hearts, etc. are used in wearable medical electronic instruments
  3. Wearable medical electronic equipment such as ECG
- two. Do not use the following persons, otherwise it may cause accidents or skin or body problems.
1. Heart disease, severe high blood pressure or
  2. Pregnant women, women who have just given birth
  3. Do not use this product in places with high humidity such as bathrooms, otherwise it may cause contacts
  4. Please read the instruction manual carefully before use, and use this product in a normal way
  5. This product cannot be consigned by air, and it is not allowed to take this product on an airplane