


Photonic Rejuvenation Beauty Instrument - DS-8811A



USER MANUAL

Remark: The device isn't waterproof, please wipe it with a wet facial tissue or wet cloth after use.

Parts Name

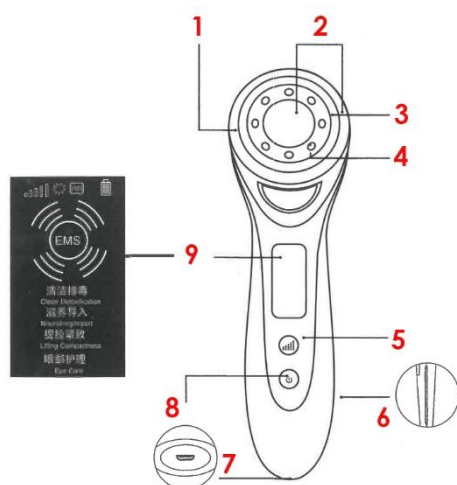

1. Makeup cotton fixing ring; 2. Massage head; 3. Lampshade; 4. LED lights; 5. Temperature adjustment; 6. ION induction strip; 7. USB; 8. Power/Mode Switch Button; 9. Operation display screen.

Functions Mode

Mode 1: Cleansing and Detoxification Function: Blue light + Vibrating Strong/Weak Gear: Vibration Strength: 1 to 5 gears

Mode 2: Nourishing Import Function: Red light + Heating + Vibrating Strong/Weak Gear: Red light and vibration changes synchronized 1 to 5 gears

Mode 3: Lifting and Tightening Function: Yellow light + Heating + Vibrating + EMS Strong/Weak
Gear: Constant Temperature + Constant Vibrating + EMS 1 to 5 gears


Mode 4: Eye Care Function: Heating + Vibrating Strong/Weak Gear: Constant Temperature + Constant Vibrating

Four Different Modes Operation

Cleansing and Detoxification(Need to use with makeup cotton)

This mode can effectively clean the dirt in the pores and make the skin bright. Use the machine to clean your face first and then use cleansing foam to clean the derived dirt. (The cleansing foam and the machine cannot be used at the same time)


1. Remove the Cotton Trap

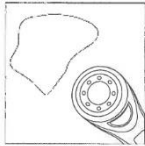

2. Put on a makeup cotton



3. Pour on the lotion

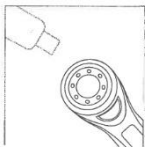

4. Long press the power button 1.5 seconds to boot the instrument



5. Hold the induction strips on both sides of the machine with your fingers



6. Facial circle cleansing (counterclockwise circle massage)



Nourishing Import(no need to use with makeup cotton)

This mode can enhance cell activity, lighten wrinkles and help skin care products moisturize the skin deeply.

1. Put a facial mask on your face or use it with skin care products.



2. Long press the power button 1.5 seconds to boot the instrument, and select the " Nourishing Import " mode.



3. Hold the instrument, touch the induction strips on both sides of the instrument, take a counterclockwise circle massage on the face to absorb.



4. Wash the remaining essence of the facial mask with water or massage to absorb according to your skin needs.



Lifting And Tightening(no need to use with makeup cotton)

This mode wakes up and relaxes the skin, helps the skin restore its radiance and vitality.

1. Apply daily skin care products or firming lotion evenly to skin.



2. Long press the power button 1.5 seconds to boot the instrument, and select the "Lifting Compactness" mode.



3. Hold the instrument, touch the induction strips on both sides of the instrument, and lifting the face from bottom to top.



Eye Care (can be used with eye cream)

This mode is mainly used for eye massage.

1. Apply eye cream evenly around the eyes



2. Long press the power button 1.5 seconds to boot the instrument, and select the "Eye Care" mode.



3. Hold the instrument, touch the induction strips on both sides of the instrument, and massage the eyes from the inside to outside.



Notice:

1. Long press the power button 1.5 seconds to boot the device, cannot be used while charging.
2. When using the device, use your fingers to hold the induction strips on both sides of the instrument to form a power circuit with the human body.
3. Use makeup remover or cleaning lotion when exporting, Don't use with cleansing foam or makeup remover when exporting.
4. The device uses micro-current ion technology, there is a slight current sense is normal phenomenon, please feel free to use.
5. The device isn't waterproof, please wipe it with a wet facial tissue or wet cloth after use.