


Color light beauty instrument

Instruction manual



- Please read the manual carefully before use.

Product structure diagram:

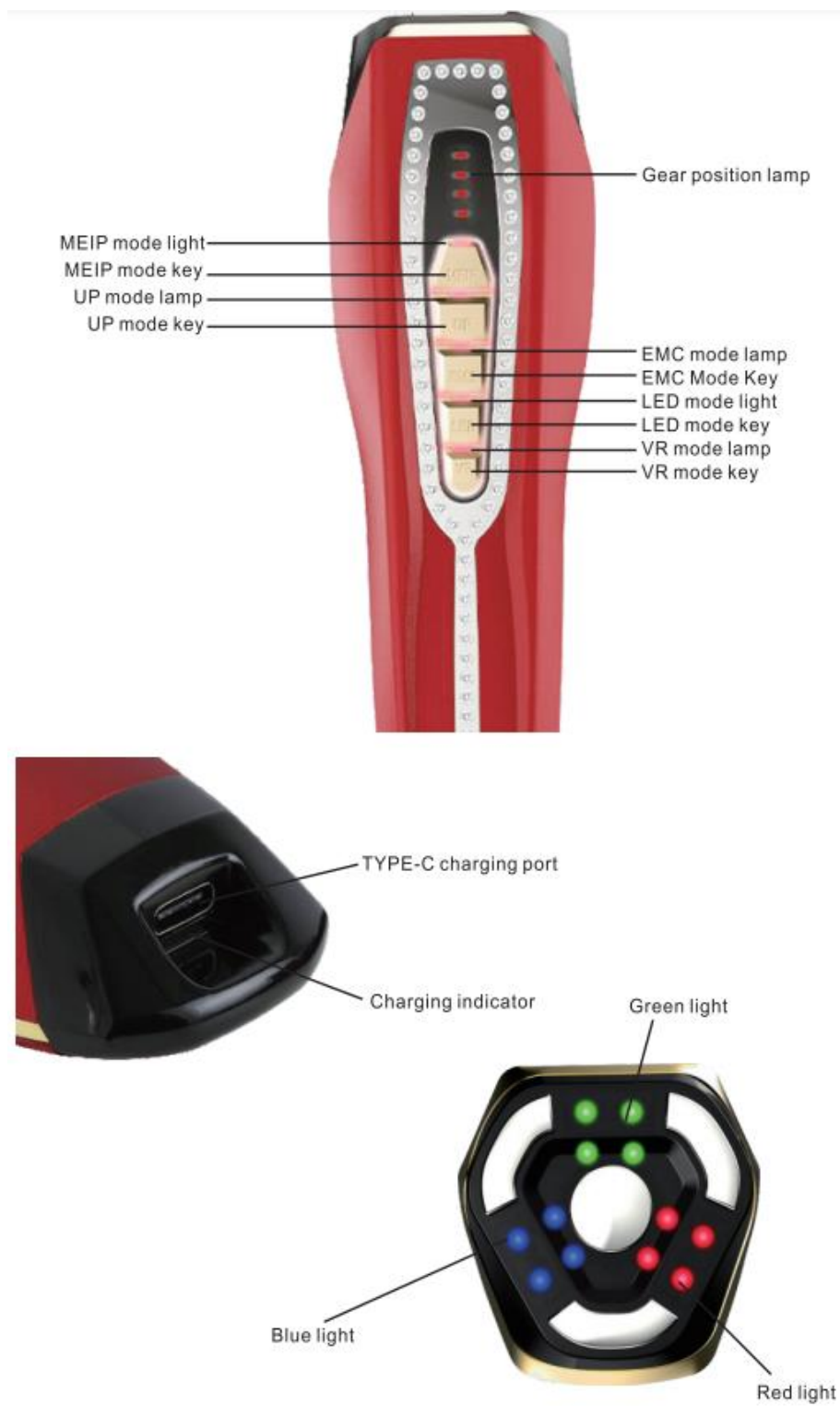


Preparations before use.

Operating instructions:

1. The machine has three modes: MFIP mode, Up mode and EMS mode; Two auxiliary modes: LED mode and VR mode.
2. Press any key of MFIP, UP, EMS and VR mode keys to start the machine and enter the corresponding mode. The auxiliary mode can be combined with the functional mode for better effect.
3. Gear selection: MFIP mode, Up mode and EMS mode are divided into four gears from weak to strong. Press the selected mode key and continue to press this key for a short time to adjust the gear. The default gear is the first gear when starting up. Press it once to increase the gear, and then press it again when it reaches the fourth gear to shut down. When it is used for the first time, whether in MFIP mode, UP mode or EMS mode, it starts in the first gear.

After starting to use, increase the strength to a comfortable gear for nursing. When reducing the strength and switching modes, press the switch of the mode used, and then adjust the gear or select the mode.

Note: After the mode is switched, the gear will automatically return to the first gear.

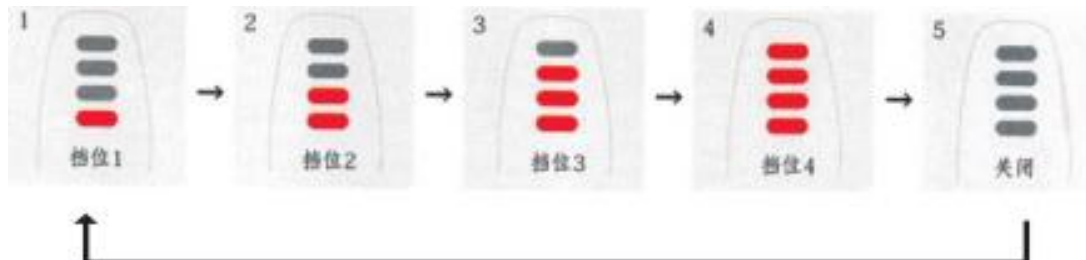

4. Select the LED mode switch, and press this key in turn. There are five phototherapy modes: blue light constant light mode, red light constant light mode, green light constant light mode, red green blue three color flashing mode and red green blue constant light mode.






Note: When the light is on, do not look at our directly look at the light source; Close your eyes when caring for the area around them.

5. Select VR vibration massage mode, and press this key to start the intermittent vibration function. The vibration massage mode can be used together with the superimposed function mode, or it can be used by opening the vibration massage mode separately.
6. Before using this product, please fully charge the product. The red indicator at the charging port is always on during charging, and the red indicator is off after charging.
7. During charging, it is forbidden to start the machine.

Pattern	MFIP mode	UP mode	EMS mode	LED mode	VR mode
Gear	<ul style="list-style-type: none"> Press this mode key to enter the first gear Press again to enter the first gear Press again to enter the third gear Press again to enter the fourth gear Press Shut Down again 	<ul style="list-style-type: none"> Press this mode key to enter the first gear Press again to enter the second gear Press again to enter third gear Press again to enter fourth gear Press again to shut down 	<ul style="list-style-type: none"> Press this mode key to enter the first gear Press again to enter second gear Press again to enter third gear Press again to enter fourth gear Press again to shut down 	<ul style="list-style-type: none"> Red green blue constant light mode The blue light is always on The red light is always on The green light is always on Bright red, green and blue 	Intermittence vibration
Function	Intermediate frequency intermittent pulse	Superpulse electroporation	Microcurrent	Light switch	Vibration massage
Mode superposition	VR mode: • The introduction of nutrient essence enables skin care products to be deeply introduced into the bottom layer of the skin, enabling the skin to deeply absorb the nutrient essence. • Use with essence, moisturizing water, facial mask and other skin care products.	VR mode: • Deeply clean the skin, remove pore dirt, and make the skin white from inside to outside. • Cooperate with facial cleanser, facial cleanser and other skin care products.	VR mode: • The effect of soothing massage and lifting and tightening can help to lighten fine lines, tighten loose skin and improve the contour. • Use with skin care products such as firming essence or firming box.	It can be used with MFIP mode UP mode, EMS mode, VR mode, etc.	It can be superposed with MFIP mode, UP mode, EMS mode, LED mode, etc.
Describe	The super pulse current can be released intermittently, and the effect can be more obvious when used first in nursing.	Simultaneous release of incoming pulse/EMS/medium and high frequency wave.	Release special pulse that can induce muscle movement.	1. Blue light mode is applicable to acne skin with strong sebum secretion and large pores, and has a good inhibition effect. The red light mode is applicable to fade fine lines and brighten the skin. Green light is applicable to improve the skin roughness, blackheads, acne, etc. caused by overwork, and effectively relieve the skin glossiness. Color light mode is applicable to people with various skin beauty needs.	It can relieve skin fatigue, and can be used in combination with other modes to have a more obvious effect and reduce current irritation.
Note: MIP mode, U mode and EMS mode cannot be mutually registered during operation. UP mode and VR mode can be overlapped with other modes. The instrument will automatically shut down when each mode is started for 10 minutes.					

Function description:

MFIP intermediate frequency intermittent pulse	Intermittently release super pulse current. The first use in nursing can make the effect more obvious.
Up super pulse electroporation	Simultaneously release the import pulse/EMS/medium high frequency wave
EMS	Release special pulse that can induce muscle movement
VR vibration mode	Equipped with vibration system to slow down current simulation.

For the first time Full light LED 	Use red, green and blue LED light energy. It is suitable for people with complex skin troubles or for anti-aging maintenance.
The second time Red LED 620nm(±10nm) 	It is suitable for people who are tired due to aging. It can recuperate bright and elastic skin.
Third time Blue LED 470nm(±10nm) 	It is suitable for people with acne and other skin problems. Create healthy and beautiful pure skin.
The fourth time Green LED 515nm(±10nm) 	Suitable for people who feel dull or uneven skin tone. Create a perfect skin with transparent impression.
Fifth time Flashing LED 	The last procedure for nursing is to help the beauty solution penetrate and absorb.

Nursing methods:



Application method of beauty solution. For full face care, please divide the face into four blocks according to the left figure, and then apply it to the whole face once and then care at the same time.

Before care:

1. Please remove the necklace, ring or any other accessories as well as all metal object on your body.
2. Please remove cosmetics, sebum, dirt, etc. In the general makeup removal and cleansing procedure before taking care.
3. According to the use mode, apply beauty liquid or condensation on the care part.

MFIP mode, UP mode --> beauty solution.

EMS mode --> massager.

- MFIP mode, UP mode = 10 minutes in total

1. Apply the beauty solution on the care part and then carry out the care.

- Do not stay at the same part for more than three seconds continuously to avoid skin problems.
- The single mode can also be used for 10 minute care.

- EMS mode = 10 minutes

1. Please apply the condensation on the care part before treatment.

- If u feel insufficient condensation during the nursing process, please add additional. Too little may cause irritation to the skin, and the nursing effect cannot be fully exerted.
- During nursing, please fully fit the four guide heads to the skin.
- Excessive care will cause skin burden. Please avoid using it more than once in a day.
- When using, if the probe surface is dirty and feels rough, please use a soft cloth or facial tissue to clean the probe surface.
- Please avoid using it when you are ill.


After care:

1. After care, gently massage the skin to help the essence remaining in the skin to penetrate and absorb. After care, please wipe the machine with a wet cloth or facial tissue to remove dirt, and then use dry facial tissue to remove moisture.

Nursing sequence:


- Please use it in a clean state after bathing.
- Please apply the beauty solution or condensation on the care part.

Upper arm: Special care, 10 minutes for unilateral arm daily care.

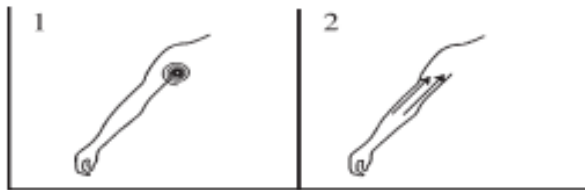

1. Draw a circle at the joint of the arm and shoulder for 2-3 circles.
2. Slide in a straight line from the elbow joint through the upper arm to the shoulder joint.

Waist: Special care, 10 minutes daily care.

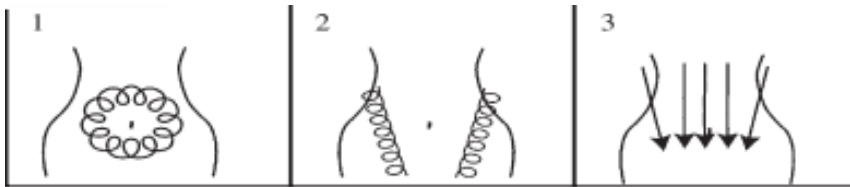

1. Slide 1 circle around the belly button.
2. Slide from the side of the waist to the groin lymph like a curve of the waist.
3. Flow smoothly from top to bottom.


Hips: Special care, 10 minutes daily care.



1. Slide upward from the joint of the hip line and thigh to the lower part of waist in a circle way.
2. Sawtooth sliding at the junction of the hip line and large clothes.

3. From the junction of the hip line and large expansion, lift and slide to the lower part of the waist in a straight line.

Expansion part (front and back): Special care, 10 minutes for unilateral leg daily care.

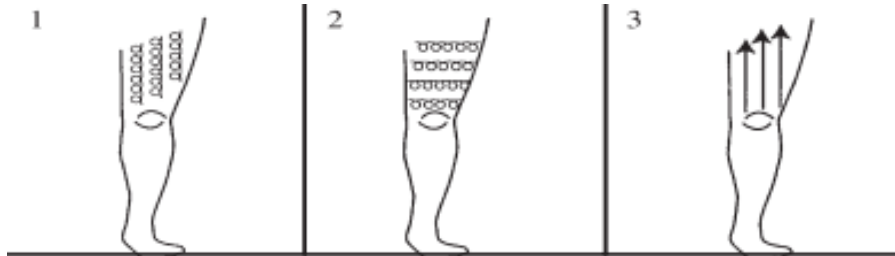


1. From above the knee (back of the knee), slide from bottom to top in a circle way to the joint of the thigh (bone line).

2. Slide laterally from above the knee (back of the knee) in a circle way, and protect it in sections from bottom to top.

3 From the top of the knee (back of the knee) to the bottom in a straight line.

- The nursing shall be performed with the same action before and after.

Calf belly: Special care, 10 minutes unilaeral leg Daily care.

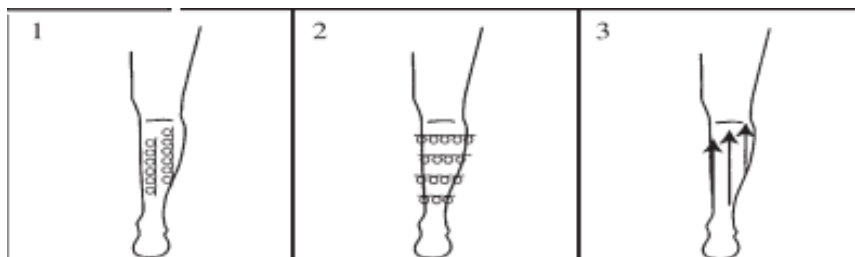

1. From the ankle to the back of the knee, slide from bottom to top in a circle to the thigh joint.

2. From the ankle to the back of the knee, slide laterally in a circle, and take care form bottom to top in order.

3. From the ankle to the back of the knee, slide from bottom to top in a straight line.

- The nursing shall be performed with the same action before and after.

Face:



1. Slide from the chin along the contour of the face to behind the ears. Then slide down to the shoulder through the auricular gland.
2. Slide from the tip of the chin along the contour of the face to the temple. Then it slides from the nose to the temple through the cheek. Finally, it slides from the head of the eye to the temple along the eye.
3. The eyebrows are lifted to the hairline and nursed from bottom to top.
 - Steps 1 to 3 shall be implemented 2 to 3 times respectively.
 - Please keep moving slowly to avoid staying at the same position.
 - If you feel that the amount of beauty solution is insufficient, please add it separately.



- Steps 1 to 5 shall be implemented 2 to 3 times respectively.
- As shown in the figure, use the instrument to care the skin of each part in sequence for about 3 seconds.
- The icon nursing procedure is the recommended nursing procedure.
- If the amount of condensation is insufficient, please add it separately.

Warning: Please avoid using it on mucous membrane, upper eyelid, eyeball, ear and around the laryngeal node.

For proper and safe use:

Be careful:

To prevent accidents or physical discomfort.

- Those who use the following medical electronic machines should avoid the possibility of accidents or physical discomfort caused by the misoperation of this electronic instrument.

Prohibit:

- Implantable medical electronic machines such as cardiac pacemakers.
- Life sustaining electronic machines such as artificial heart and lung.
- ECG scanner or drip mounted medical electronic machines.

Note: This product is a household electronic product and cannot replace medical devices or drugs.

Warning:

Do not use it for the following persons:

- When suffering from heart disease, malignant tumor or abnormal blood pressure.
- Suffering from diseases, infectious diseases and acute diseases accompanied by fever.
- During dental treatment.
- Patients with immune deficiency diseases such as AIDS. Who have decreased immunity or used immunosuppressants.
- Hemophilia disease, hemorrhagic disease, diabetes, facial neuralgia.
- Those with abnormal perception and pain.
- Crab foot swelling, constriction, facial melanosis, light allergy, abnormal pigment.
- Long term use of steroids or micro blood tension caused by liver dysfunction.
- Allergic constitution, atopic dermatitis, desmatosis, sensitive skin or skin allergy.
- Have a history of epilepsy or feel unwell due to flashing light.
- People with sequelae in the head or eyes due to traffic accidents.
- Those who have undergone eye surgery (those who have undergone laser surgery should consult a doctor before proceeding).
- During pregnancy, physiological period and lactation.
- Young children or those unable to express their own will.
- People with alcohol and hangover.
- People taking drugs that induce drowsiness.
- People who feel extremely tired or in poor physical condition.
- People with skin problems such as spots and inflammation caused by chemicals.
- Dermatitis caused by or ever caused by cosmetics.
- People who have recently received chemical skin rejuvenation treatment in skin care salons.
- People who have skin problems after using beauty instruments.









Do not use in the following parts:

- Outside the designated position.





- Pericardiac, eyeball, upper eyelid, lips or middles mucous membrane of mouth.
- Location of pathological stain, nevus or viral wart.
- Parts causing skin tingling and inflammation due to excessive sunlight.
- Parts with pathological stain due to itching, fever, physical stimulation, etc.
- The position, pudenda, anus, around the laryngeal knot, and ears that have undergone surgery.
- Parts subject to hair removal or shaving recently.

Note: Patients who are being treated by doctors should consult with them before use.

Precautions during use:

 Prohibit	<ul style="list-style-type: none"> • Please do not use it together with other beauty instruments. • Please do not use it for any purpose other than beauty, and do not use it when it is damaged. This beauty instrument is a special beauty care device. Please do not use it for treatment or other purposes. • Do not look at or directly look at the LED light. • Do not use it on the same part continuously, so as not to damage the skin or cause skin problems such as infection, inflammation, bleeding, pigment precipitation or worsen the symptoms. • Do not use on pets.
 Prohibited use of water field	<ul style="list-style-type: none"> • Do not use or keep it in the bathroom or in a place with heavy humidity. Do not place it in a damp place with water (on the washstand, etc.).
 No wet water	<ul style="list-style-type: none"> • Do not soak, spray or pour liquid on the instrument. (May cause electric shock, fire or instrument failure)
 No wet hands	<ul style="list-style-type: none"> • Do not touch the instrument with damp hands.
 Prohibit	<ul style="list-style-type: none"> • Please do not use it near flammable objects such as banana water volatile oil polishing water spray. • Do not use charging wires other than those specially used.
 Unplug	<ul style="list-style-type: none"> • When unplugging the charging cable, be sure to hold the front end plug and then unplug it. Do not pull the wire to remove the plug.
 Prohibit	<ul style="list-style-type: none"> • Do not damage the charging cable. Damage, forced bending, pulling, twisting, placing heavy objects above, clamping wiring, etc. (May cause fire, electric shock or instrument failure)
 Must guard	<ul style="list-style-type: none"> • Only applicable to household sockets • When inserting the plug, please firmly insert the charging wire.

WARNING:





 Unplug	<ul style="list-style-type: none"> • The machine body or charging line emits smoke, abnormal heat, odor and deformation
 Prohibit	<ul style="list-style-type: none"> • The instrument may be hot during use, which is normal, and there is no doubt about safety. • Do not connect this instrument to metal or electrical products. (Mobile phone, TV, etc.)
 Decomposition prohibition	<ul style="list-style-type: none"> • Do not disassemble, reconstruct or repair the instrument by yourself. If you need to repair the instrument, please contact the purchase shop.
 Must guard	<ul style="list-style-type: none"> • Before using this instrument, please be sure to remove metal ornaments, glasses or contact lenses. Metal jewelry: earrings, necklaces, watches, etc • For first-time users or those not adapted, please set the power as the most primary use. • Please close your eyes before using when caring the skin around your eyes. • If you feel abnormal skin or body pain during use, please immediately stop using and consult a doctor.

WARNING:

Precautions for use of this instrument, please read the following items carefully before using.

1	Special care for each part=once a week ●Take care twice a week in the first two weeks. The third week after that depends on the state of the skin, about once a week. ●Excessive care may have adverse effects on the skin. Daily daily care is carried out at the frequency of every part=once a day or twice a day ●Please use it according to the skin condition.
2	●Facial care should be carried out after removing makeup, cleansing, removing makeup and sebum dirt, and applying special beauty solution.
3	The skin condition may vary depending on the physical condition. ●Please pay attention to the skin condition during use. If the skin has abnormal symptoms such as rash, redness and itching, please stop using and consult a specialist.

WARNING:

 Unplug	<ul style="list-style-type: none"> • Before cleaning the instrument body, please unplug the charging line from the instrument and socket.
 Prohibit	<ul style="list-style-type: none"> • Do not install metal parts or dust them on the charging line or plug. • Do not drop or hit the instrument.
 Prohibit	<ul style="list-style-type: none"> • Please keep or use it in a place where children cannot get it. • Please do not store it in the following places. • Places with direct sunlight, heavy humidity and dust dispersion • Places with extreme temperature, such as high temperature or below zero, and vehicles exposed to high temperature and sunshine • Places where there may be a risk of ignition or exposure to wind and rain • Places where chemicals, gases or violent vibration occur • Near heating equipment or heat source, or around magnetic substances
 Must guard	<ul style="list-style-type: none"> • Cleaning method • Please use a soft cloth towel soaked in water and wrung dry. If a coarse cloth towel is used for wiping, the instrument may be scratched or the function may be weakened. ⌘ Do not use alkaline solvent, alcohol, banana water, volatile oil and other solvents, or polishing powder, hot water, or chemicals to clean the machine. ⌘ Do not use the instrument when it is stained with special water.

Please check before use:

- If the instrument has not been used for more than half a year, be sure to check whether the instrument is abnormal before use.

Disposal of rechargeable batteries:

1. The built-in rechargeable lithium battery of this product contains substances that may pollute the environment. When discarding this product, please do not put it together with general domestic waste, and take out the battery.
2. Before disassembling the product and taking out the lithium battery, be sure to turn off the power or cut off the power.
3. Take out the lithium battery. Push the snap on the upper cover of the product, open the upper cover of the product, and then take out the battery.
4. The lithium battery taken out should be safely disposed and put into the battery recycling center designated by the government, which is conducive to environmental protection.

Special note: Take out the battery only when discarding this product. Make sure that the battery has been exhausted before taking it out. Non professionals are not allowed to remove it.