


EYE MASSAGE INSTRUMENT



An instruction manual

Product mix

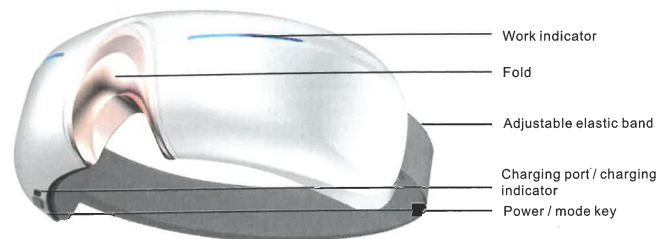

Product parameters

Product Name	Eye massage instrument	Rated voltage	5V
Charging mode	TYPE-C	Rated power	5W
Thermal temperature	42-45°C	Control mode	One key control
Executive standard	GB4706.1-2005 GB4706.10-2008		

Matters needing attention

Please read this manual carefully and operate according to the instructions.

1. It can not be used when there are symptoms of eye disease. It is prohibited for patients with major eye diseases such as eye surgery, glaucoma, cataract and retinal detachment.
2. It is forbidden to wash the product with water. Please use cotton soaked with alcohol for cleaning.
3. Do not disassemble this product by yourself.
4. This product has a hot surface. People who are not sensitive to heat must pay attention to it when using it.
5. Minors, the disabled and the elderly shall be used under the supervision of adults.
6. This product is a personal health massage product and cannot replace drugs and medical devices.
7. It is prohibited to disassemble, impact, squeeze the internal lithium battery or put it into the fire; In case of severe swelling, please do not continue to use.
8. Do not use the battery after soaking in water.
9. There is a lithium battery inside, please do not put it in the high temperature environment.
10. For sensitive skin, facial allergic dermatitis and plastic materials on the face, use this product under the guidance of a doctor.
11. If you feel unwell during use, please stop using immediately.

Operating instructions

Precautions before use

1. Please remove the contact lenses, hair clips, hair ornaments, etc. before use.
2. Please use it with your eyes closed, put the product on your eyes and then turn it on to avoid damage to the airbag due to overpressure.
3. Do not use at high temperature and humidity.
4. After the use, the eyes are blurred. After closing the eyes for 1-3 minutes, the eyes will gradually become clear.

Charging method

Insert the type-C data cable into the charging port of the product. The red indicator is always on when charging, and the green indicator is always on after charging. The product cannot be started during charging.

Usage time

The default working time of this product is 15 minutes. If you need to extend the service time, you need to start it again.

Usage method

This product is divided into 4 types: vibration massage ds-b10, vibration hot compress ds-b11, luxury voice ds-b12, supreme Bluetooth ds-b15

How to use vibration massage

1. Long press the power / mode key to start the machine, and open the free mode. The green indicator light at the power key is always on.
2. Short press the power / mode key to switch gears, and the third gear will cycle in sequence.
1st gear: start vibration
Gear 2: start vibration
Gear 3: vibration reinforcement
3. Press the power / mode key for a long time to shut down the machine, and the indicator light goes out.

Application method of vibration hot compress

1. Long press the power / mode key to start the machine, and open the free mode. The green indicator at the power key
The light is always on, and the blue indicator light of the fuselage is on.
2. Short press the power / mode key to switch gears, and the 4th gear will cycle in sequence.
Gear 1: start vibration and hot compress
Gear 2: start vibration and hot compress
Gear 3: vibration reinforcement and hot compress
Gear 4: turn on hot compress and turn off vibration
3. Press the power / mode key for a long time to shut down the machine, and the indicator light goes out.

How to use the luxury voice model

1. Long press the power / mode key to start the machine, and the voice prompts "welcome to use the eye massager". The green indicator light at the power key is always on, and the blue indicator light of the body is on.

2. Press the power / mode key for a short time to switch the mode and give a voice reminder. The default free mode is startup. The machine is cycled and switched in free mode, eye protection mode, vibration mode, decompression mode and sleep mode.

Free mode: vibration + music + air bag + hot compress

Eye protection mode: Music + air bag + hot compress

Vibration mode: vibration + music + airbag

Decompression mode: Music + airbag

Sleep mode: air bag + hot compress

3. Press the power / mode key for a long time to shut down the machine. The voice prompts "thank you for using" and the indicator light goes out.

How to use supreme Bluetooth

1. Long press the power / mode key to start the machine, and the voice prompts "welcome to use the eye massager". The green indicator light at the power key is always on, and the blue indicator light of the body is on.

2. Press the power / mode key for a short time to switch the mode and give a voice reminder. The default free mode is startup. The machine is cycled and switched in free mode, eye protection mode, vibration mode, decompression mode and sleep mode.

Free mode: vibration + music + air bag + hot compress

Eye protection mode: Music + air bag + hot compress

Vibration mode: vibration + music + airbag

Decompression mode: Music + airbag

Sleep mode: air bag + hot compress

3. Bluetooth Music: you can also use the mobile phone Bluetooth to connect with the machine, and then operate the music player of the mobile phone to play music on the machine.

4. Press the power / mode key for a long time to shut down the machine. The voice prompts "thank you for using" and the indicator light goes out.

Daily maintenance

- When the surface of the massager is dirty, wipe it with a clean soft cloth dipped in a little water. Do not wipe with volatile liquid such as alcohol. The product cannot be disassembled and cannot be rinsed in water.
- If it is not used for a long time, please cover the product to avoid dust.
- Water stop or other liquids shall be prevented from entering the machine during cleaning to avoid short circuit during use.

Troubleshooting

Serial number	Problems encountered	Possible causes	Processing method
1	The sound of "wheezing" and "dada" appears in the action	The normal sound generated by air pump inflation, motor rotation, air valve deflation cloth sleeve friction, etc	No processing required
2	Automatic shutdown after startup	Insufficient lithium power supply	Please charge
3	Unable to turn on power	Insufficient lithium power supply	Please charge
		Charging	Please disconnect the power supply for use
4	The indicator does not light up when charging	Poor contact of USB socket	Please unplug the USB cable
		Poor power contact	Ensure good contact
5	The air pressure is too weak	Insufficient lithium power supply	Please charge

If the above processing is still abnormal, please turn off the power and contact the dealer you purchased.

Do not disassemble and repair the product without authorization. The company will not be responsible for product damage and personal injury.

Disposal method of rechargeable battery

1. The built-in rechargeable lithium battery of this product contains substances that may pollute the environment. When discarding this product, please do not put it together with general domestic garbage and take out the battery.
2. Before disassembling the product and taking out the lithium battery, be sure to turn off the power or cut off the power.
3. To take out the lithium battery, pry the buckle of the upper cover of the product, open the upper cover of the product, and then take out the lithium battery.
4. The lithium batteries taken out should be safely disposed of and put into the battery recycling center designated by the government, which is conducive to environmental protection.

Special note: only take out the battery when discarding this product. Please ensure that the battery has been exhausted before taking it out. Non professionals should not remove it.

Warranty Card

Product information

Model: _____

Colour: _____

Purchase channel: _____

Customer information

Full name: _____

Address: _____

Purchase date: _____

Certificate

Este produto é inspecionado para atender aos padrões nacionais de qualidade e é entregue

Product Name: Eye massage instrument

Product model: _____

Product color: _____

Quality testing: _____

Date of manufacture: _____

Product executive standard: _____

