


Smart Fat Removal Machine – 9150 / 9151



User's manual

IMPORTANT NOTE:

This product is available in two versions, the 9150 is the simplest (only pulse and vibration) and the 9151, which additionally has a hot compress function; please read the instructions that correspond to the model you have purchased. thank you!

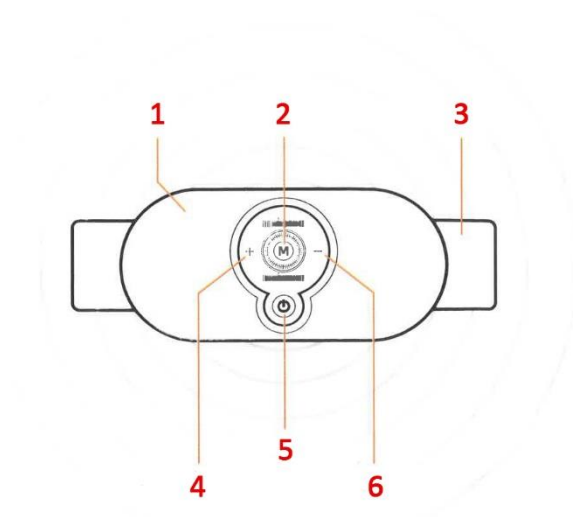

1.- Host; 2.- MODE key; 3.- Belt; 4.- Boost key; 5.- POWER button; 6.- Weaken key.

Product specifications

Product model: 9150

Product Name: Smart Fat Removal Machine

Rated power: 3W

Charging time: about 2 hours

Charging method: USB (DC 5V1A)

Rated voltage: DC 3.7V

Rated time: 15min

Gear mode: 9 levels of vibration intensity, 18 levels of pulse intensity Features: dual-purpose waist and abdomen, vibration massage + pulse massage

Executive standard: GB4706.1-2005 GB4706.10-2008

Product model: 9151

Product Name: Smart Fat Removal Machine

Rated power: 5W

Charging time: about 2 hours

Charging method: USB (DC 5V1A)

Rated voltage: DC 3.7V

Rated time: 15min

Gear mode: 9 levels of vibration intensity, 18 levels of pulse intensity Features: dual-purpose waist and abdomen, vibration massage + pulse massage, warm compress

Executive standard: GB4706.1-2005 GB4706.10-2008

Operation method

9150 operating instructions:

Open the packing box and take out the fat-removing machine, first charge it with the charging cable that comes with it, and then unplug the charging cable before it is fully charged.

- Long press the "POWER" key to turn on the instrument, and starts to work.
- Press "MODE" to switch between 5 pulse modes: 1. Mixed mode 2. Low-frequency beating mode 3. High-frequency beating mode 4. Micro-point pulse therapy mode 5. Acupuncture mode.

- Short press the "POWER" key, the vibration intensity of the motor will cycle from 1-9 levels in turn.
- After selecting the massage mode, press the "+" key to increase the intensity, and the "-" key to decrease the intensity, which has the function of increasing or decreasing the intensity. There are 5 pulse modes with 18 levels of intensity.
- During use, long press the "POWER" key to turn off the instrument.

9151 operating instructions:

Open the packing box and take out the fat-removing machine, first charge it with the charging cable that comes with it, and then unplug the charging cable before it is fully charged.

- Long press the "POWER" key to turn on the instrument, and starts to work.
- Press "MODE" to switch between 5 pulse modes: 1. Mixed mode 2. Low-frequency beating mode 3. High-frequency beating mode 4. Micro-point pulse therapy mode 5. Acupuncture mode.
- Long press the "MODE" heating key to turn on/off,
- Short press the "POWER" key, the vibration intensity of the motor will cycle from 1-9 levels in turn.
- After selecting the massage mode, press the "+" key to increase the intensity, and the "-" key to decrease the intensity, which can increase or decrease the intensity. There are 5 pulse modes and 18 levels of pulse strength.
- During use, long press the "POWER" key to turn off the instrument.

How to connect power

Please use the USB charging cable that comes with the instrument to connect the mobile phone charger or power bank (The output voltage is 5V, the current is not more than 1.5A), The red light will be on after power-on, and it will go out when it is fully charged.

Why doesn't the product have a charging adapter

The product supports universal Android USB data cable Support the standard universal mobile phone adapter and power bank on the market. (Note: output voltage is 5V, output current> 500mA, recommended output current> =1A)

Safety Precautions

The purpose of the warning signs and illustrations in the manual is to enable you to use the product safely and correctly, and to prevent harm to you and others. The warning signs, legends and their meanings are as follows

Product features and functional characteristics

1. Designed according to ergonomics, beautiful and scientific, in line with the neck curve
2. One-key smart low-frequency electric pulse:

3. Using 3D intelligent fitting technology, adjust the massage according to the curve of the abdomen:
4. The comfort of the mixed massage function is further improved.

Waste battery disposal

Before the product is discarded, the battery must be removed from the product. When removing the battery, the product must be powered off. First use a screwdriver to unscrew the shell. After removing the shell, the battery can be removed. After removal, the waste battery should be sent to the battery recycling department for disposal or put into the battery recycling bin.