

# INSTRUCTION MANUAL FOR THE MUSCLE STIMULATOR

## DMAF0184



Please keep and read it (safety Attention) before using the caring items hereby is recorded necessary important contents for safety This notice divides into three classes when you make a mistake with handling and using.

### **Drawing indication:**

**Warning:** contents that include the possibility of fatality

**Caution:** Contents that include the possibility of bodily injury or property damage.

# DANGER

## Prohibition:

Prohibition to people who are using medical treatment like below:

- The pulse adjustment machine is for medical treatment by electronic and etc that inset into the body.
- The artificial heart is for medicinal treatment by electronic and etc that maintains living.
- An electrocardiograph is for medical treatment by electronic machina and etc that are by paste.

## What is a low-frequency Treatment?

Our body occurs to the lowest electricity, which we cannot feel. but it is getting sensitiveness action by outside electronic impetus, and according to it, which take several changes on body. By body character, let body take a minimum electricity on special part of body, to exercise muscular by impetus for relieving fatigue and dissolution stress.

- Relieving symptoms for tightening shoulder, waist and back
- Betterment of muscular paralysis
- Massage effect
- Betterment of diet

## Product

### Product features

1. The low frequency treatment is by rhythm and soft impetus.

The low frequency is a special electricity with regular rhythm or stop. According to suitable power and speed on body by low frequency, which make a rhythm and sort impetus on body. The softness of low frequency are as same as massage effect like pressure, pinch and percuss by human hand.


2. Choice by 6 modes.

1) Choice 6 modes with simple operation

2) It will get off power automatically after 12 minutes from stop, when you forget power off.

# How to use

## How to operate

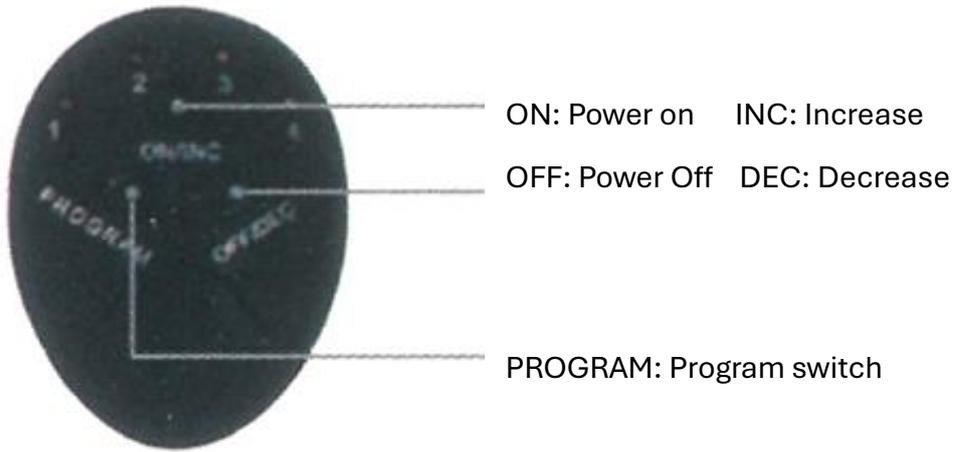

## Lamp for mode indication



When it presses the "power on" and "rising intensity" button and the lamp comes to dawn little by little.

When it presses the "power off" and "down intensity" button and the lamp comes to dark little by little

- Button for both "power on" and "rising intensity", "power off and "down intensity" intensity divides with 10classes
- Button for mode choice: choice by 6 modes (A—F)
- Dry battery cover: which is to rear open the cover and change the dry cell battery.

## How to use the Pad:

- 1) Always wash the area of your body where you intend to use the adhesive gel pads.
- 2) Pull the pads off the vinyl cover by gripping the pads

3) Attach the other adhesive gel pad to a desired part of the body(arms, shoulders, etc.)

1. Don't attach it with "power on" because it could show you.

2. Don't touch attach side of adhesive pad, it make decreasing using time.

### **Procces for operation:**

1. Take out the equipment from the case.  
Please start it with at lowest output when it is beginning.
2. Place the gel pads at desired areas of the body.
3. Please push switch for power on.
4. Choice mode.
5. Adjustment for output level.
6. Working time is 12 minutes per one.
7. To push switch for "OFF/DEC" button until let lamp off. when it doesn't act during working time.

### **How to use an effective**

1. The low frequency is suitable to 12 minutes one time and 2-3 times a day, infused it a long time or feel harder with impetus, then you will feel tired at moment.
2. It's able to avoid 30minutes after dinner, and favor to do after shower or warm washing.

### **Taking care after using**

1. Please turn off power after using.
2. Please separate adhesive gel pads from machine.
3. Please attach adhesive gel pad to pad holder.  
(To prevent adhesive pad for drying, which must paste it pad holder)
4. The dirty point etc. on cine should be cleaned with wet cloth by neutral cleaning material and wiped it by drying cloth again.
5. Don't keep in sunshine, hot temperature and humidity, wet and dust places, closer fire, or shocking places etc.
6. To keep it away from children who cannot touch.

7. Don't touch drying battery terminal with wet hands.
8. If you don't use it for a long time. (more than 1 month) please separate battery from machine and store it at another place.

### **Main machine**

1. The dirty points etc. on the case are cleaned with wet cloth by neutral cleaning material and wipe it with drying cloth again.
2. Do not keep it in high or low temperatures or humidity when storing the unit.

### **Adhesive gel pad**

1. You will feel uncomfortable if the pasting side of adhesive pad gets dirty or pasting power is getting weakness more.

Then, you can use wet towel and clean dirty side for regaining pasting again at moment, but, if it gets moisture too much, that will re-act to pasting.

2. The adhesive gel pads(electrodes) should be replaced with the same type and model.

### **Change the battery**

1. Please open the cover.
2. Please confirm "+" and "-" and insert correctly.
3. Please close the cover.
4. Please turn on switch and confirm indication lamp

## Install instructions



1. Open the battery shell



2. Install the battery(2\*AAA batteries) correctly



3. Assemble the host device on abdominal pad



4. Peel off the adhesive pads.



5. Assemble the host device on arm pads



6. Attach the pads on your arms.